# Certificates

- 1. Medical Electrical Equipment: General Requirements for Safety IEC/EN 60601-1.
- 2. Top level standard for electromagnetic compatibility for electrical medical equipment IEC/EN 60601-1-2.
- 3. Requirement for safety of hospital bed such as construction, mechanical strength, mechanical hazards, protection against fluid ingress and essential performance IEC/EN 60601-1-2-52.

# Features

	ES-12HD	ES-12DW	ES-12HDM	ES-12HF
Size	L:2220 * W:1000 * H:450 mm			
X-ray	Ο	Ο	Option	Option
CPR and Quick release	Ο	Ο	Option	Option
Trendelenburg	0	0	Option	Option
Angle indicator	Ο	Ο	Option	Option
Central brake and steering system	Ο	Ο	Option	Option
Weight scale	X	Ο	X	X
Manual Back-up	X	X	Ο	X

# JosonCare

Factory:

No.11, Aly. 77, Ln. 180, Hushan St., Guishan Township, Taoyuan County 33348, Taiwan(R.O.C.) Tel:886-3-3290925 Fax:886-3-3290921

Tel:886-3-3290925 Fax:886-3-329092

Office

No.280, Chengdu Rd., Xitun Dist., Taichung City 40751, Taiwan(R.O.C.) Tel:886-4-23160743 Fax:886-4-23160740

> Http://www.hospitalbed-josoncare.com/ E-mail:overseas\_ma@joson-care.com.tw

> > Document No. : JCQB-078 Release Date :2013.07.01 Version : A

# JosonCare

# Electric Weight Scale Bed ES-12DW



Specialty • Safety • Simple • Reliable





# ES-12DW

# ICU Electric Hospital Bed

**ES-12DW** is made by Joson-care excellent technology and follow the medical standard as a design idea. The design not only enhances the patients' comfort and safety, but also improves the convenience for caregivers when operating. Moreover, the innovative design ideas and the patents can meet all the requirement during the nursing care process.

#### Safety And Firm construction

The ES-12DW ICU Hospital Bed's frame is made of steel with powder coating and the mattress' platform is made of ABS plastic which has ventilation holes.






## X ray


There is an X-ray cassette under the back board allowing examinations to be performed directly on the bed without the risks of patient mobilization.

### Drop-down Side-rails

Both sides have the patient control panel which can manipulate Hi-Lo, back, knee and Auto-chair independently and the outside panel is used by caregivers to prevent the spread of infection.



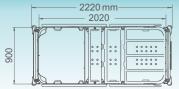

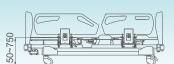





#### Over Bed Length: 2220mm

Over Bed Width(siderail up): 1000mm

Over Bed Width(siderail stored):
 950mm

Over Mattress Base Length: 2020mm

# Specifications

Over Mattress Base Width: 900mm

 Over Bed Height Adjustment : 450~750mm

Back rest: 0°~75°

Knee gatch: 0°~38°

Trend/Reverse Trend tilt: ±12°

# 3 Angle Indicator



Nursing staff can monitor the angle of Trendeleburg & Rev.

Trend from the tilt indicator and monitor the angle of back-rest.



## 4 Quick release

Release the back-rest for emergency



## 5 Mattress Stopper



# The New Design of Head/

## 6 Foot Board's Lock

There are two locks located at outside of head/foot board which

are much easier for removing.





## Attendant Control Panel

#### Weight scale

1. Weight capacity: 0~250kg

2. Minimum measurement unit: 100g

3. Kg/Lb switch mode

4. Patient monitor mode

5. Hold mode

Control panel which can manipulate Hi-Lo, back, and knee

The control panel can be locked to prevent unexpected contact



Hold function — Patient leave alarm — Weight variation alarm — Prevent patient falling alarm —

Electrical Trendelenburg and Rev. Trendelenburg

Electrical CPR
Release all positions and height to be flat position.